### SQL Error Code 20 - Database not Attached

Last Modified on 10/30/2025 2:00 pm PDT

#### **Summary**

This document shows you how to fix a SQL Error Code 20 for a database that is not attached.


#### **Before you Begin**

- 1. Make sure the Windows account you are logged into has Administrator privileges
  - · Check Administrator Privileges if you are unsure.
- 2. IMPORTANT: Make sure Windows OS is fully updated before continuing.
  - To check if you are fully updated, go to the Microsoft support article here:
    https://support.microsoft.com/en-us/windows/install-windows-updates-3c5ae7fc-9fb6-9af1-1984-b5e0412c556a

# How to Attach your SQL Database

### 1. Check to see if Database File is in Proper Location

- Open Windows File Explorer
- Navigate to C:\ProgramFiles\Microsoft SQL Server\MSSQL14.PBI\_MYOVISION\MSSQL\DATA
- Confirm MyoVision and MyoVision\_log database files are in the folder

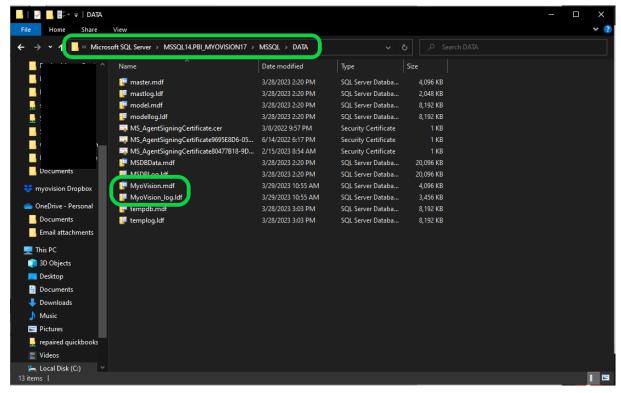

### 2. Move, Attach, or Troubleshoot database files

• If the Database file is in the DATA folder;

- Move to **Troubleshooting** section
- If the Database file is not in the DATA folder,
  - Perform Step 9 and 10 Moving and Attaching Database Files in the Install SQL Server 2017 support article
  - If a permission error appears, like the one shown below, follow the steps under <u>Troubleshooting</u>




# **Troubleshooting**


If you are getting a **permissions issue error**, you will need to change the permissions for the DATA folder. Once you've done this, go back to the previous section and try steps 1 and 2 again.

### 1. Open DATA Folder Properties

 Open File Explorer and navigate to C:\Program Files\Microsoft SQL Server\MSSQL14.PBI\_MYOVISION\MSSQL

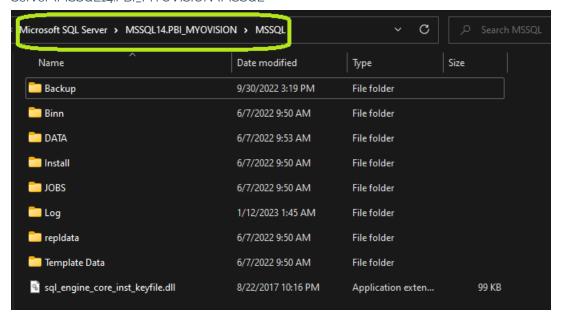


• Right-click DATA, select Properties

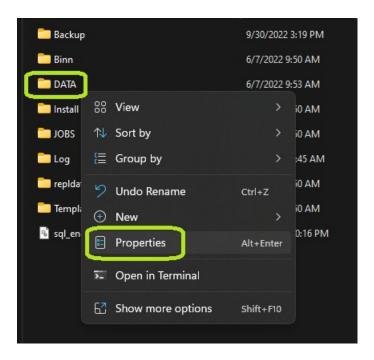


## 2. Edit Groups/User Names


• Click the Security tab in the window that pops up.

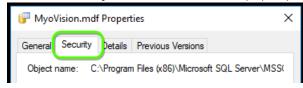


- Check under "Group or user names" to see if Users is listed.
  - If so, skip ahead to Step 3: Edit Permissions
  - If Users is not listed, click **Edit**

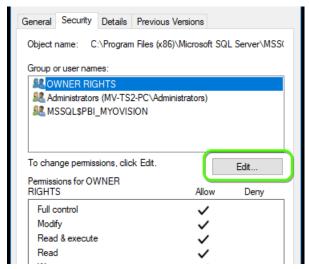


• In the Permissions window that opens, click Add

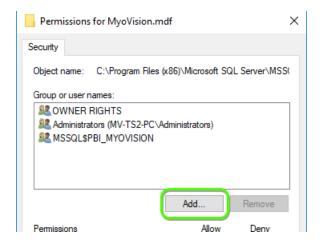


• In the Select Users or Groups window, click **Advanced** 

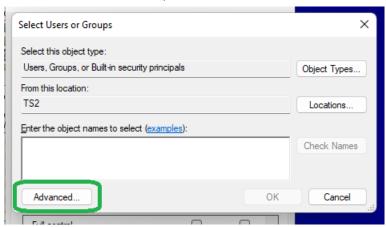


- Click Find Now, find and selected the "User" account in the list at the bottom of the window
  - If there is no User account, look for Domain User or other group that ends in "User"
- Click Ok, then click OK in Select Users or Groups window

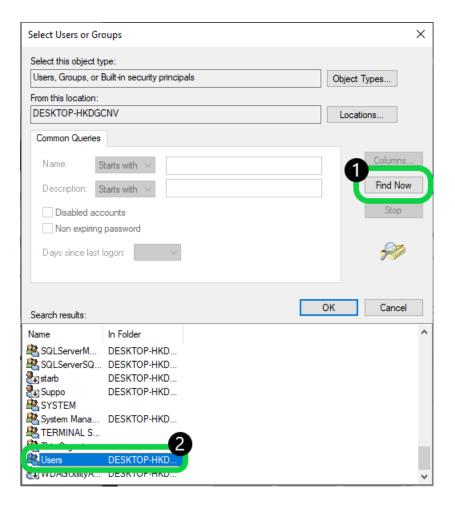

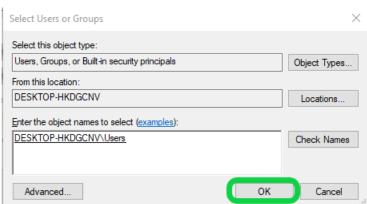

## 3. Edit Permissions

• Now that Users or Domain Users has been added, click on it under "Group or user names" and allow **Full Control**. Click OK to finish

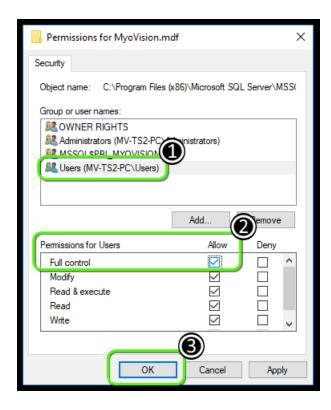

### 4. Attach the database

- Perform Step 10. Attaching Database Files in the Install SQL Server 2017 support article
  - If the database is not working, an SQL error will appear. See the links below for help with the SQL Error
    - SQL Error 10
    - SQL Error timeout
    - SQL Error 'sa' login failed

**Note:** If you are still having trouble getting your MyoVision database to work properly, please feel free to contact our technical support at 1-800-969-6961, option 2 (8:30am - 4:00pm, Monday - Friday) (+1 206-448-3464 International)